



THE No. 1 System  
for Solvent Safety

# Safety Funnel with hinged lid "Marco"

Smart Protection when collecting Liquid Waste!



# Clever Technology, Strong Protection!

Safety Funnel "Marco"

✓ Improved Closing

Comfortable handling: easy, quick opening and closing with "click"-closure.

✓ Hinged Lid with Sealing Ring

For secure closure, when the funnel is not in use.

✓ Stable Hinge

A double-screw hinge, constructed of rust-free stainless steel: for long life and stability.

✓ Optimized Dirt Sieve

For collecting mixing magnets or large dirt particles. Laboratory bottles can also be placed here to dry off. The hole size has been optimized, in order to improve the running-off of chemicals.

**The sieve can be removed for cleaning.**

✓ Various Screw-Threads 

Suitable for different sizes of collecting canister (see Table on Page 3).

✓ Safety Lance 

Ensures for safe and good flowing of sensitive chemicals, without any "stuttering".

## Dimensions

\*Installation height:  
+ 105 mm

Lid open vertical:  
+ 270 mm




# Configurations

Safety Funnel "Marco"



Fig.	Screw-Thread of container	Lance	Electrically conductive (incl. grounding cable and clamp)	Integrated Level Control (with Float in mm)	Delivery with Adapter
<b>A</b>	GL 45	-	-	-	-
-	GL 45	-	•	-	-
<b>B</b>	GL 45	•	-	-	-
-	S 50	-	-	-	•
-	S 50	•	•	-	•
-	S 50	•	-	-	•
-	S 51	-	-	-	•
-	S 51	•	•	-	•
-	S 51	•	-	-	•
-	S 55	-	-	-	•
-	S 55	-	•	-	•
-	S 55	•	•	-	•
-	S 55	•	-	-	•
-	S 60/61	-	-	-	-
<b>C</b>	S 60/61	-	•	-	-
<b>D</b>	S 60/61	•	•	-	-
-	S 60/61	•	-	-	-
<b>E</b>	S 60/61	-	•	• 150	-
-	S 60/61	-	•	• 190	-
-	S 65	-	-	-	•
-	S 65	-	•	-	•
-	S 65	•	•	-	•
-	S 65	•	-	-	•

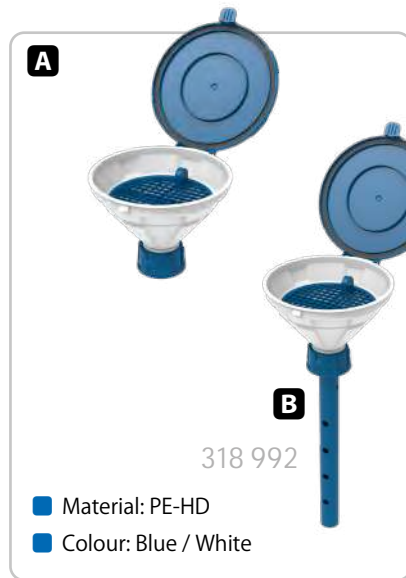


**Information about "Delivery with Adapter":**  
 The safety funnel "Marco" is available in various container screw-thread sizes. In order that we may produce in an environmentally-friendly and economic manner, certain variants are delivered together with a screw-thread adapter. Those to which this applies can be seen in the column "Delivery with Adapter".








## S.C.A.T. Europe GmbH

Waldecker Straße 7  
64546 Mörfelden-Walldorf  
Deutschland

 +49 (0) 61 05 / 30 55 86 - 0


 +49 (0) 61 05 / 30 55 86 - 99


 [info@scat-europe.com](mailto:info@scat-europe.com)

 [www.scat-europe.com](http://www.scat-europe.com)

## S.C.A.T. Americas, Inc.

1704 East Blvd., Suite 101  
Charlotte, NC 28203  
USA

 +1.704.438.9477

 +1.704.373.2603

 [info@scat-americas.com](mailto:info@scat-americas.com)

 [www.scat-americas.com](http://www.scat-americas.com)

 **Via Torre, 2**  
 **28040 Borgo Ticino (NO) - Italy**  
 **+39 0321 036125**  
 **[customercare@vivaxsrl.it](mailto:customercare@vivaxsrl.it)**  
 **<https://www.vivaxlab.com>**

